


# cosine<sup>®</sup>

艾柯医疗  
ACCUMEDICAL


## Distal Access Catheter

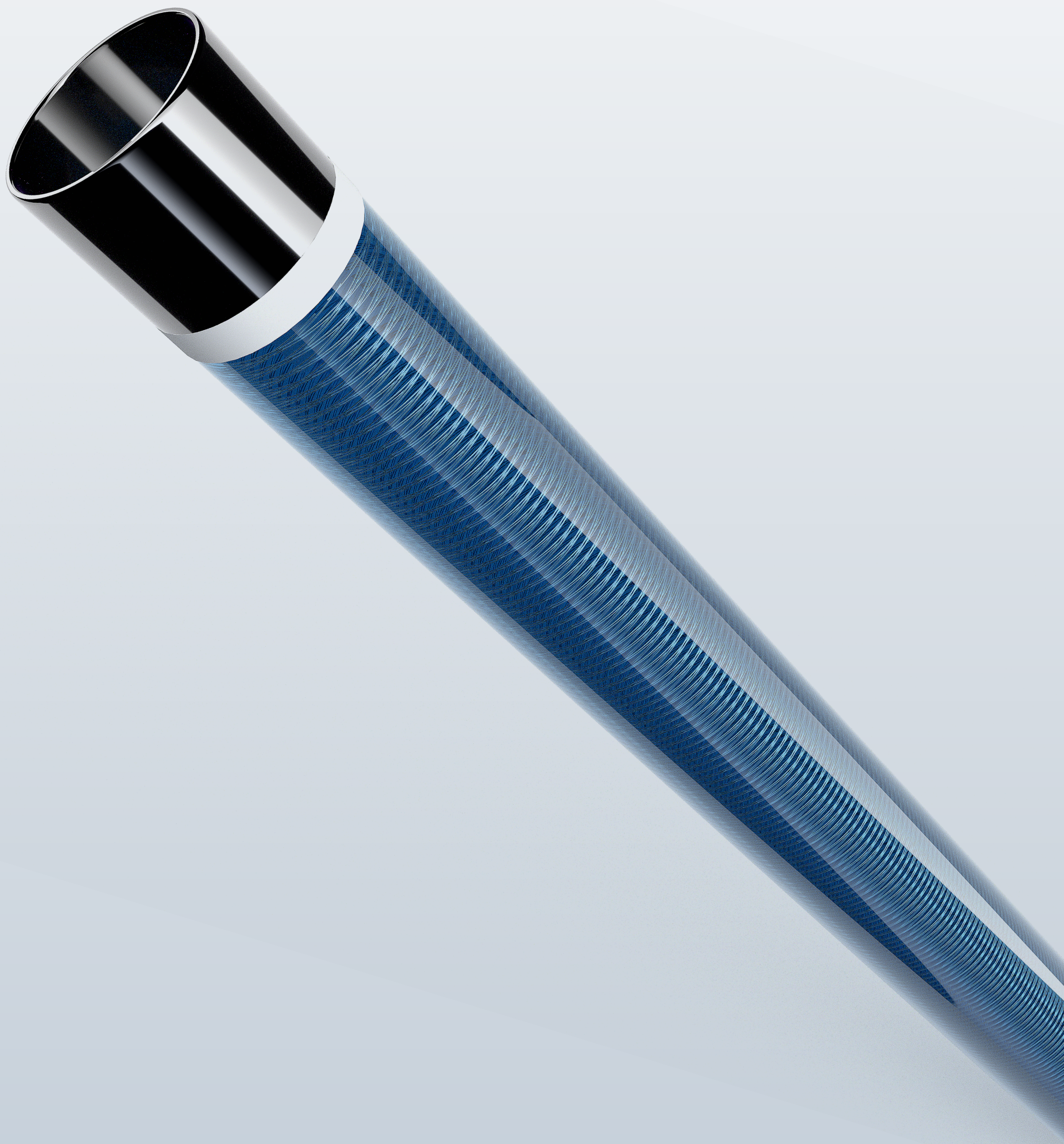

Smooth, steady and compatible

Medical Device Registration Certificate No. of P.R. China: 20223030188

# | Strong Support

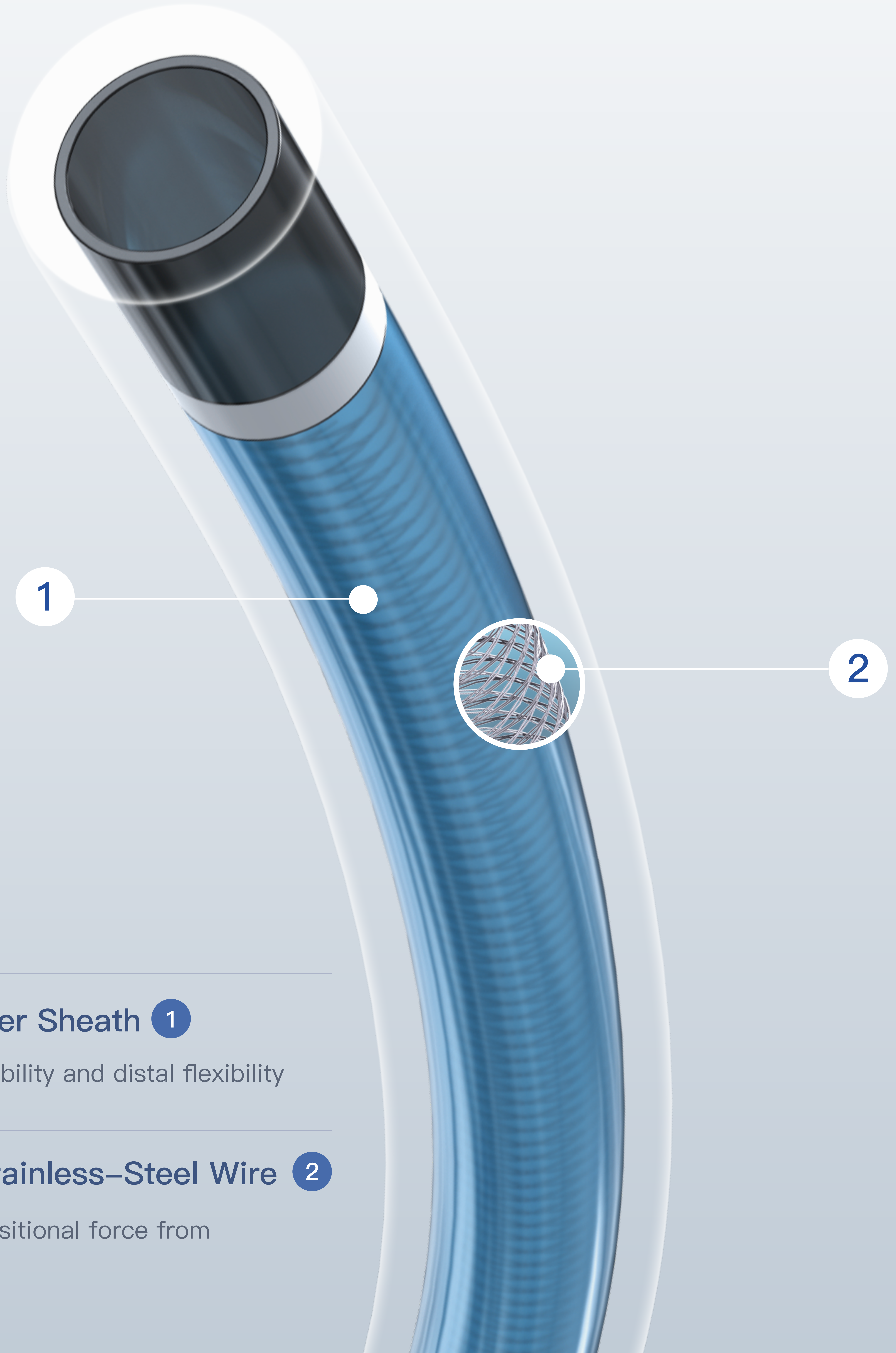

---

## 10-segment Outer Sheath 1


Provides proximal stability and distal flexibility

---

## Twin-Braided Stainless-Steel Wire 2

Provides smooth transitional force from proximal to distal

# | High Compatibility

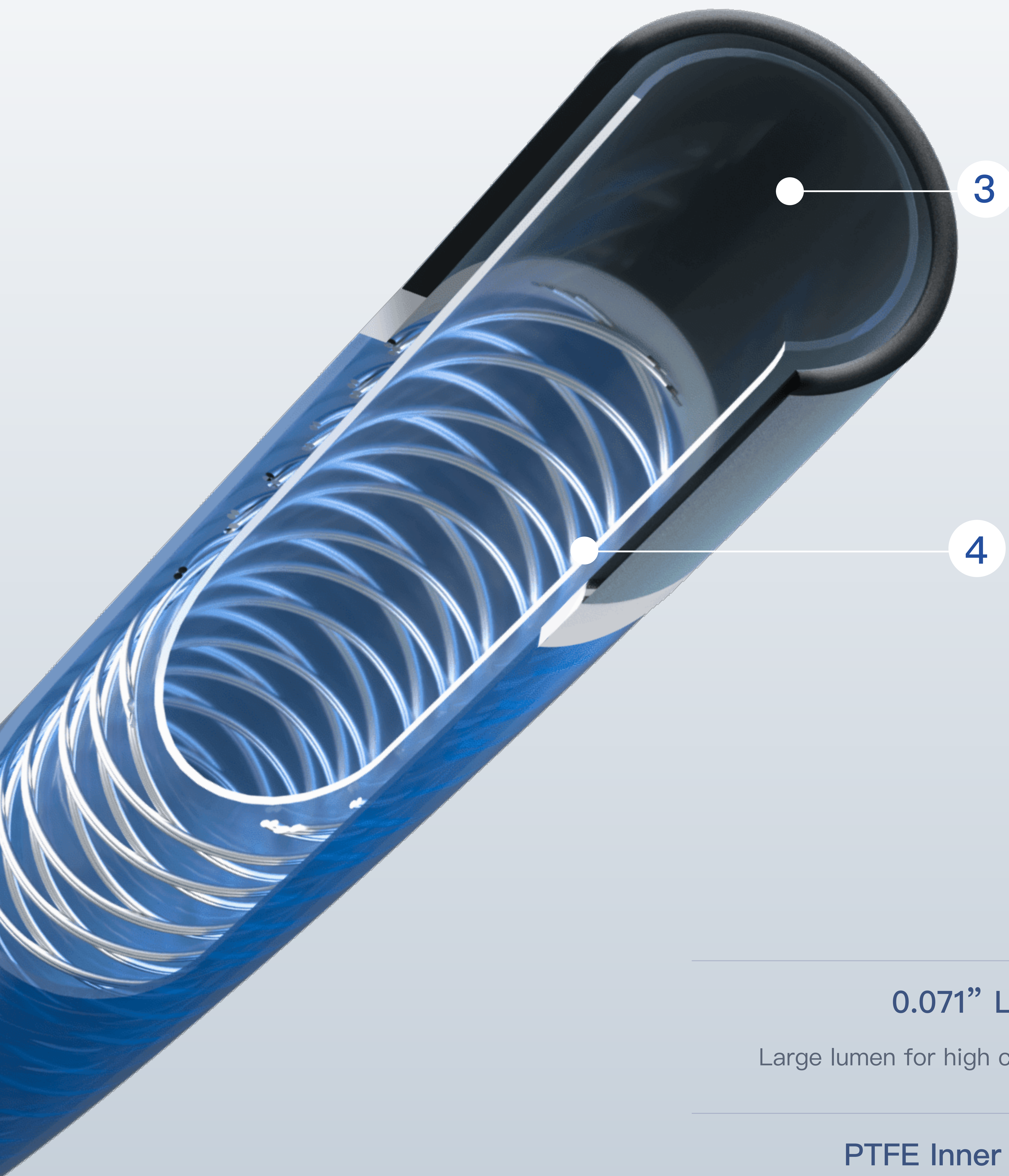

---

0.071" Lumen 3

Large lumen for high compatibility

---

PTFE Inner Layer 4

Smooth delivery of microcatheters and  
interventional devices

# | Distal Access

---

## Distal Tip with Radiopacity 5

3mm macro-molecule distal tip, less stimulation to the vessels than metals

---

## Distal End Soft Zone 6

18–21cm length, bringing support point closer to treatment site, and increasing stability for device delivery and deployment

---



## 40cm Hydrophilic Distal End Coating 7

Reduce resistance



cosine<sup>®</sup> Distal Access Catheter

Specifications

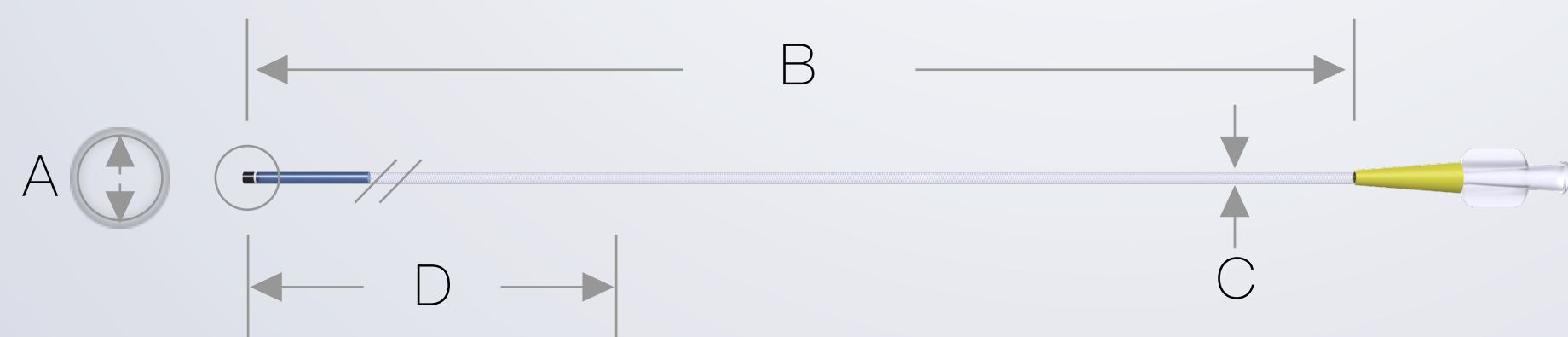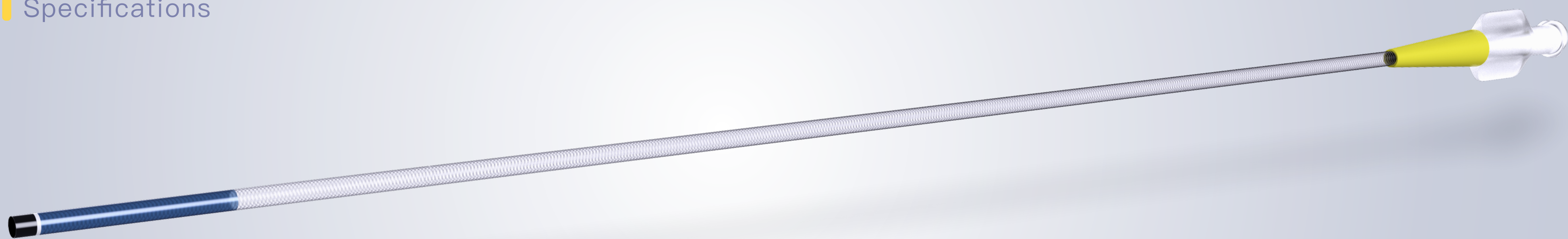

A: Inner Diameter (ID): 0.055" (5F, 1.4mm) / 0.071" (6F, 1.8mm)  
B: Length: 70–138 cm  
C: Outer Diameter (OD): 0.068" (1.73mm) / 0.084" (2.13mm)  
D: Hydrophilic Coating Length: 40 cm

Medical Device Registration Certificate No. of P.R. China: 20223030188

| SKU           |          | OD     |        | ID     |       | Effective Length |
|---------------|----------|--------|--------|--------|-------|------------------|
| Support<br>5F | SC502SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 75cm             |
|               | SC504SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 85cm             |
|               | SC506SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 95cm             |
|               | SC508SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 105cm            |
|               | SC510SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 115cm            |
|               | SC516SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 125cm            |
|               | SC521SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 130cm            |
|               | SC526SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 135cm            |
|               | SC529SHR | 0.068" | 1.73mm | 0.055" | 1.4mm | 138cm            |
| Support<br>6F | SC601SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 70cm             |
|               | SC603SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 80cm             |
|               | SC605SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 90cm             |
|               | SC608SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 105cm            |
|               | SC610SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 115cm            |
|               | SC611SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 120cm            |
|               | SC616SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 125cm            |
|               | SC621SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 130cm            |
|               | SC626SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 135cm            |
|               | SC629SHR | 0.084" | 2.13mm | 0.071" | 1.8mm | 138cm            |
| Soft<br>5F    | SC501SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 70cm             |
|               | SC503SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 80cm             |
|               | SC505SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 90cm             |
|               | SC506SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 95cm             |
|               | SC508SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 105cm            |
|               | SC510SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 115cm            |
|               | SC511SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 120cm            |
|               | SC516SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 125cm            |
|               | SC521SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 130cm            |
|               | SC526SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 135cm            |
|               | SC529SHX | 0.068" | 1.73mm | 0.055" | 1.4mm | 138cm            |
| Soft<br>6F    | SC601SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 70cm             |
|               | SC603SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 80cm             |
|               | SC605SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 90cm             |
|               | SC607SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 100cm            |
|               | SC610SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 115cm            |
|               | SC616SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 125cm            |
|               | SC626SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 135cm            |
|               | SC629SHX | 0.084" | 2.13mm | 0.071" | 1.8mm | 138cm            |

AccuMedical Beijing Inc.  
1210, Building A Pacific Century Place, 2A Gong Ti Bei Lu, Chaoyang  
District, Beijing, CHINA



www.accumedical.com  
+86 010 6506 5018